

Ultrapure water for the daily needs.

Easy and economical! TKA MicroLab.



**ASTM Type II water.
Instantly and fresh.**

TKA WATER
PURIFICATION
SYSTEMS

TKA MicroLab. All that you need.



■ **ASTM Type II Ultrapure water.**

For example, for

- Rinsing laboratory glassware
- Preparing and diluting buffers, reagents, tissue culture media and dyes
- Sample preparation for analytical methods such as flame AAS

■ **Instantly and fresh.**

With TKA MicroLab.

This is the password for converting normal tap water to ASTM Type II ultrapure water.

All in one! TKA MicroLab is a compact complete-system that produces 6 litres of ultrapure water per hour. The built-in 6 litre tank for ultrapure water has a conical bottom for complete emptying and is equipped with a sterile vent filter and a pressure pump for supplying equipment connected downstream (e.g. autoclaves, analyzers). If required, a flexible dispenser with a sterile filter can be connected to the tank. The system has a small footprint, is very attractively designed and looks good in any laboratory.

TKA MicroLab operates with the same high-tech purification procedures as the well known TKA Pacific systems.



Outstanding design!

TKA MicroLab – Not only the elegant exterior advances laboratories into the future!

TKA MicroLab.
Ultrapure water at the press of a button! Position it, connect it, use it!



Economical!

Two long-life cartridges. Each requires replacement only when really necessary!



One, two, three...

... rapid cartridge change due to quick-connects.




Space saving!

Simply stand it on the laboratory bench, or mount it directly to the wall.



Convenient!


Swivel-mounted display - always easy to read

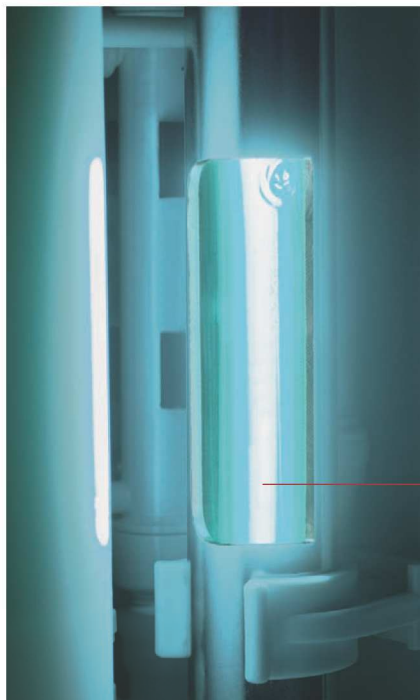

UV-lamp!

Option of a UV-lamp for micro-organism reduction.

TKA MicroLab Flow chart.

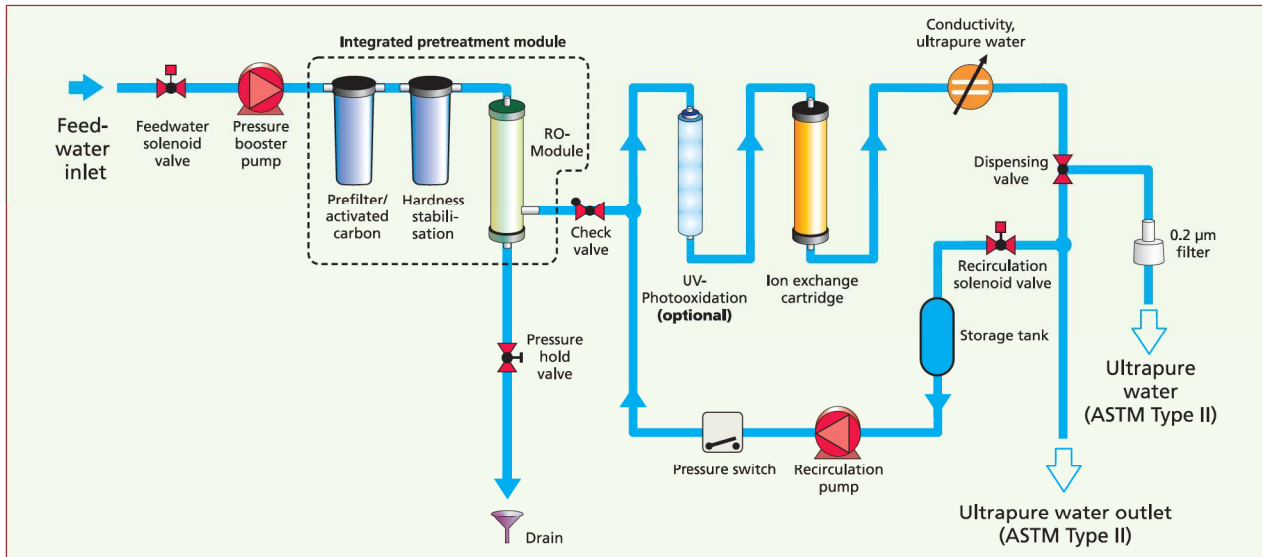

The result: Ultrapure water of the required quality:

ASTM Type II

6 Litres per hour. Always freshly produced and ready to be taken at any time - either from the built-in 6 litre tank or via the separate dispensing valve with sterile filter.

Conductivity in µS/cm:	0.067 – 0.1
Resistance in MΩ×cm at 25 °C:	15 – 10
Bacterial content in CFU/ml:	< 1
Particle content:	< 1
Permeate performance per hour:	6 litres
Flow rate per minute:	1.0 litre



Absolutely sterile!

Ultrapure water is dispensed through a standardly supplied sterile filter.

More ideas = more benefits

How TKA MicroLab differs from other compact-class systems:

Powerful performance!

- Permeate flow 6 l/h
- Flow from sterile filter outlet 1.0 l/min
- 6-Litre tank with integrated pressure and recirculation pumps

Lasting economy!

- Water purification progresses through separate modules, each with quick connects for simple individual replacement

• Module 1

Pretreatment module with integrated reverse osmosis membrane

• Module 2

Cartridge containing high-quality ultrapure resins with very long service intervals

Easiest to operate!

- Just press a button to get ultrapure water - no adjustments necessary, the microprocessor control controls itself
- Ultrapure water quality is ensured by the limiting value input
- RS 232 Interface
- The display can be tilted to your optimal reading position

Placement anywhere

- Stand it on the laboratory bench
- Mount it on the wall

Usable worldwide!

- Whether 220/230 V, 110/120 V or other voltages - and even with fluctuating voltages - the automatic voltage regulator continually supplies the system with a constant 24 volts

Almost maintenance-free!

- The tank is made of ultrapure water resistant, pigment-free and food-safe polyethylene and has a conical bottom outlet for complete emptying
- Can be easily cleaned and disinfected





Specifications

	TKA MicroLab	TKA MicroLab UV
Permeate performance at 15 °C:	6 l/h	6 l/h
Flow rate:	1.0 l/min	1.0 l/min
Operating pressure in bar, min/max:	1 – 6 bar	1 – 6 bar
Conductivity:	0.067 – 0.1 µS/cm	0.067 – 0.1 µS/cm
Resistance at 25 °C:	15 – 10 MΩ×cm	15 – 10 MΩ×cm
TOC Value:	-	< 30 ppb
Bacterial content in CFU/ml:	< 1	< 1
Particle content:	< 1	< 1
Supply voltage:	Automatic voltage regulation to 24 V	Automatic voltage regulation to 24 V
Power consumption:	0.06 kW	0.06 kW
Connector size:	R 3/4"	R 3/4"
Dimensions, W x D x H in mm:	305 x 400 x 545	305 x 400 x 545
Weight:	22 kg	22 kg
Article no.:	08.0060	08.0061*

* with UV-lamp

Consumables and accessories

Wall-mount	Article no. 09.2212
Pretreatment module with RO-membrane	Article no. 09.2006
Ultrapure water cartridge	Article no. 09.1018
Sterile filter	Article no. 09.1003
Vent filter for tank	Article no. 22.0091

Distributed by:



Carl Stuart Limited

ADVANCED APPLIED TECHNOLOGIES

Contact Us:

Irl Ph: 01 4523432

UK Ph: 08452 30 40 30

Web: www.carlstuart.com

Email: info@carlstuart.com



TKA Wasseraufbereitungssysteme GmbH

Stockland 3

D-56412 Niederelbert

Telephone: +49 (0) 26 02 / 10 69 9-0

Telefax: +49 (0) 26 02 / 10 69 9-50

eMail: info@tka.de

www.tka.de