

CertoClav Vacuum Pro Series

LABORATORY AUTOCLAVE

>> OUR
ALL-IN-ONE
AUTOCLAVE <<



- ✓ 8 – 22 liters
- ✓ Media temp. sensor
- ✓ Vacuum technology
- ✓ Cloud documentation
- ✓ Programmable
- ✓ Class B, N, S cycles



Distributed by:

Lab Unlimited
CARL STUART GROUP

Tallaght Business Park
Whitestown, Dublin 24,
Ireland
D24 RFK3

Tel: (01) 4523432
Fax: (01) 4523967
E-mail: info@labunlimited.com
Web: www.labunlimited.com

Quatro House, Frimley Road,
Camberley,
United Kingdom
GU16 7ER

Tel: 08452 30 40 30
Fax: 08452 30 50 30
E-mail: info@labunlimited.co.uk
Web: www.labunlimited.co.uk

CERTOCLAV
Since 50 years.

Catalog 2017

CertoClav Vacuum Pro Series

SUITABLE FOR ANY APPLICATION




- Microbiology, Bio-Chemistry**
Sterilize culture media precisely with the integrated media sensor
- Food industry**
Food preservation
- Material testing**
Perform hundreds of cycles in succession, fully automatically without user interaction
- Pharmacy**
Sterilize hollow instruments, porous materials, packaging, plastics, rubber and textiles
- Universities**
Robust, ideally for any application and training operation



CERTOCLAV
Since 50 years.

CertoClav Vacuum Pro Series

STAND-ALONE OR SYSTEM INTEGRATED

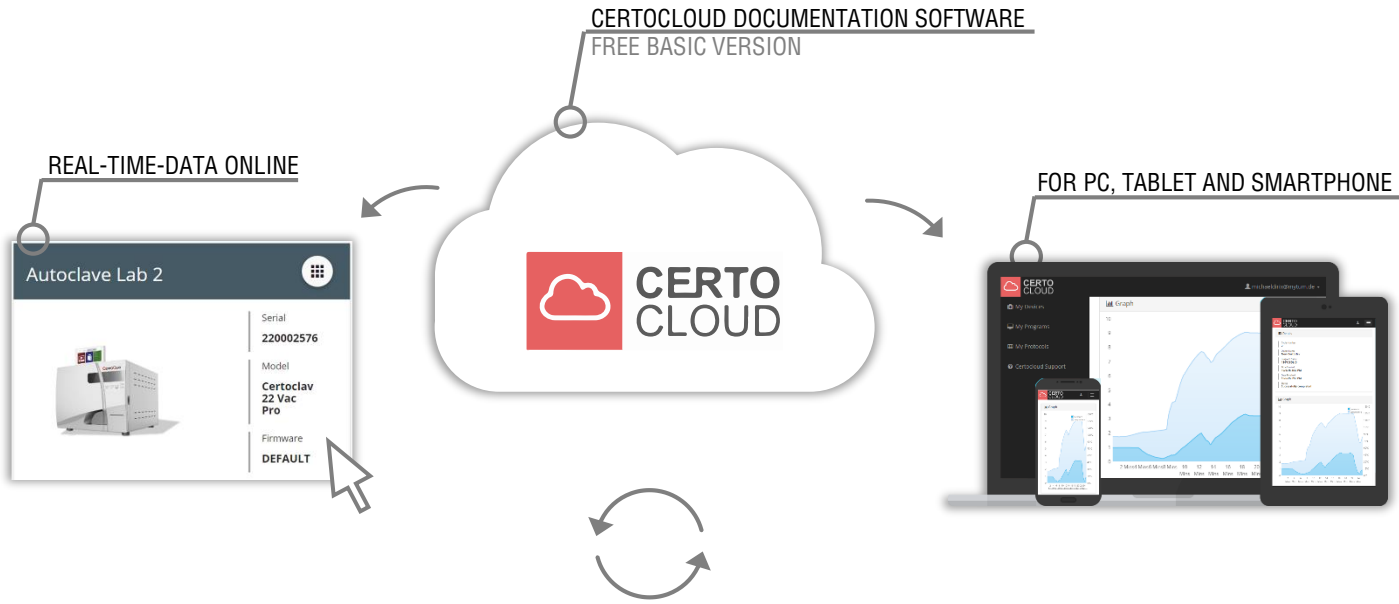

Connect the autoclave via LAN



Connect the autoclave via W-LAN



Receive Email notifications



Print protocols with the protocol printer



Label sterilized items with tracing information



Export protocols to USB flash drive



Export protocols to SD CARD.



Track the complete sterilization process via scanner.



Stand alone or fixed water connection

CertoClav Vacuum Pro Series

STAND ALONE – OR – FIXED WATER CONNECTION

AUTOCLAVE ORDER CONFIGURATIONS

Order number	Chamber size: 8 liters	Chamber size: 12 liters	Chamber size: 18 liters	Chamber size: 22 liters	Internal water tanks (Stand alone)	Fixed water connection. Including R/O- System for automatic refill
8500486	✓				✓	
8500486W	✓					✓
8500487		✓			✓	
8500487W		✓				✓
8500488			✓		✓	
8500488W			✓			✓
8500489				✓	✓	
8500489W				✓		✓

ACCESSORIES



Reverse osmosis system

Automatic water supply for the autoclave. Requires tap water connection. Power: 110-220V, 50-60Hz. Order number: 8500499



Filters for reverse osmosis system

Order numbers:
8500510 (Sediment filter, replace every 6 months)
8500511 (Precarbon filter, replace every 6 months)
8500512 (Membrane filter, replace every 30 months)



Tray set

Order numbers:
8l: 8500478
12l: 8500479
18l: 8500480
22l: 8500481



Tray for bottles

Order numbers:
8l: 8500478A
12l: 8500479A
18l: 8500480A
22l: 8500481A



Label printer

Print expiration date and protocol information on self-adhesive label and mark your goods. Order number: 8500502



Barcode scanner

Track the complete sterilization history of an item. Order number: 8500503



Service set

Includes door sealing, air filter, safety valve, quick-connect drain hose. Order number: 8500495



External printer

Print protocols with a printer. Order number: 1300187



CertoCloud premium license


- Enable SMS notifications
 - Storage for 100.000 protocols
 - Premium Remote support
 - Real-time-data online
- Order number: 8500212

CertoClav Vacuum Pro Series

INTUITIVE USER INTERFACE

TOUCH SCREEN


✓ Log into your account

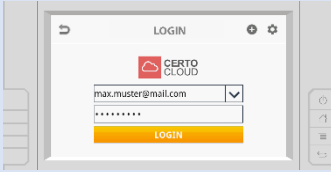


✓ Choose programs




✓ Manage your programs

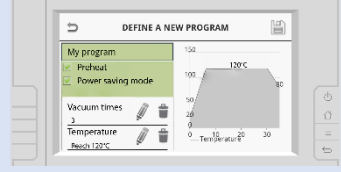

✓ Contact support

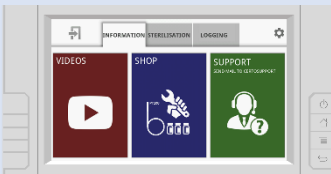

✓ Watch protocols

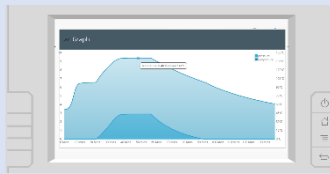


✓ Follow the integrated step by step video assistant



FREE AUTOCLAVE MANAGEMENT SOFTWARE


✓ Create an user account to manage your CertoClav products.

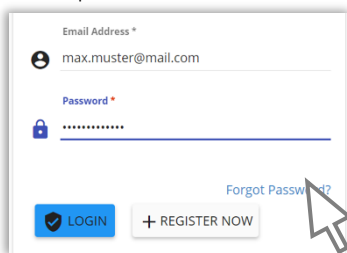


✓ Access real-time data from your autoclave from all over the world.




✓ Configure SMS and email notifications.

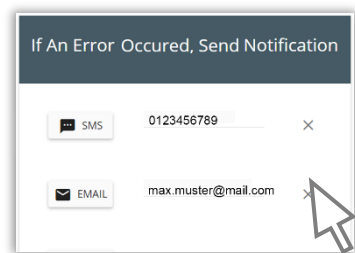


✓ Scroll through protocol history. Export protocols to signed PDF, Excel or ZIP.

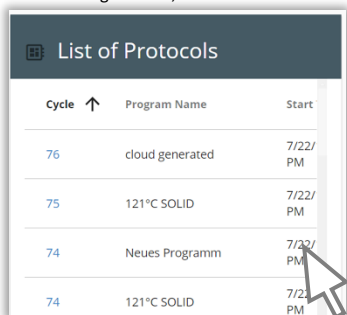

✓ Show all protocols from the autoclave graphical or as a list.

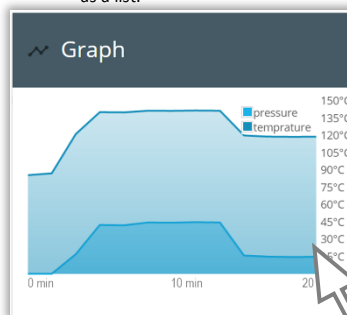


✓ Confirm the authenticity and validity of printed out protocols with a few clicks.

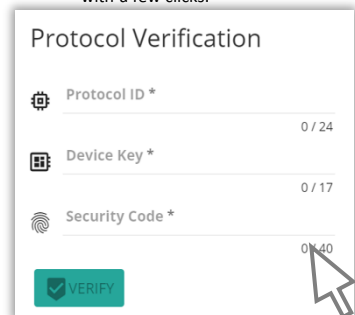

CertoClav Vacuum Pro Series


STERILIZE SOLID MATERIALS AND LIQUIDS

STERILIZATION OF SOLID MATERIALS

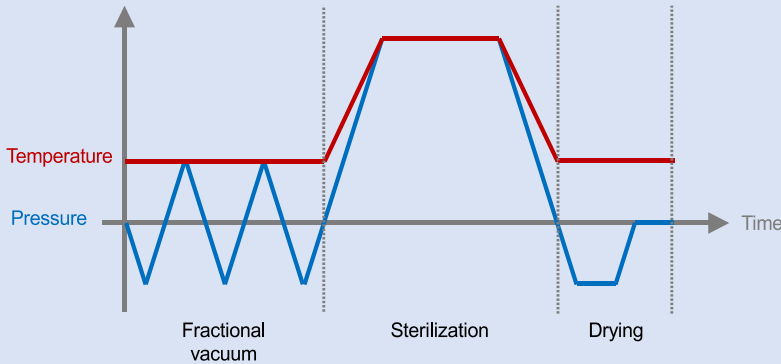

- ✓ All sterilization programs for solids follow a typical Class-B cycle.
- ✓ Precise combination of vacuum, pressure and temperature for optimal sterilization of solids.
- ✓ The fractionated vacuum presses the saturated water vapor through the sterilization package and into each microscopic gap.

STERILIZATION OF LIQUIDS WITH MEDIA SENSOR

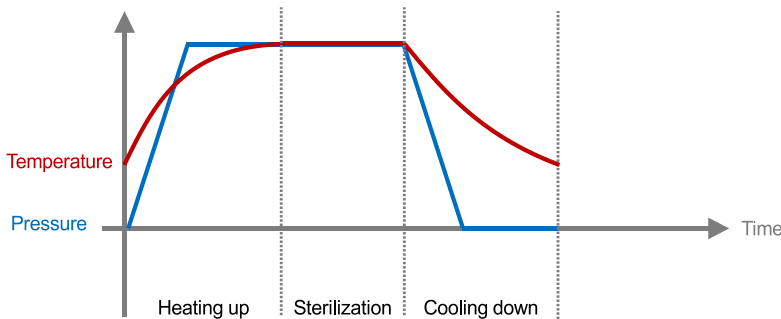

- ✓ All sterilization programs for liquids follow a heating phase with venting followed by a sterilization phase and a final cooling phase..
- ✓ Sterilize liquids at 121°C for 20 minutes.
- ✓ The integrated media sensor measures the real temperature of liquids and ensures an optimum sterilization result.

STERILIZEABLE MATERIALS



Instruments



Unpacked materials, single and double packed instruments



Liquids in bottles




Liquids with media temperature sensor



Foodstuffs



Textiles

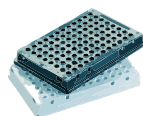

Plastics, rubber materials



Porous materials, sand



Hollow materials



Filters, waste

CertoClav Vacuum Pro Series

CUSTOM PROGRAM AND AUTOMATED MATERIAL TESTS

PROGRAMS

	Vacuum times	Temperature [°C]	Holding time [min]	Drying time [min]	Total time [min]
Solid	1	121	20	3	32
Solid quick	1	134	4	3	20
Porous	3	121	20	7	50
Porous quick	3	134	4	7	27
Hollow	3	121	20	10	47
Hollow quick	3	134	4	10	34
Liquid	0	121	20	0	45 to 120
Prion	3	134	20	10	50
B-D test	3	134	4	7	31
Vacuum test	1	-	15	-	15
Clean process	-	-	-	-	29
User-Defined	1 to 10	105 to 134	1 to 60	1 to 40	

MATERIAL TESTS

Standard Program

Run for 100 cycles automatically

START



Set repeat function for material tests



Each sterilization program can be executed as often as desired



Watch real-time-data on www.certocloud.com



























Receive a notification, after the material test has been finished



CertoClav Vacuum Pro Series

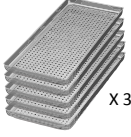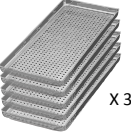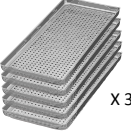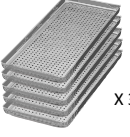STERILIZE SOLID MATERIALS AND LIQUIDS

LOADING CAPACITY

Volume of chamber	8 liter	12 liter	18 liter	22 liter
Diameter of chamber	17 cm	20 cm	25 cm	25 cm
Depth of chamber	32 cm	36 cm	35 cm	45 cm
50 ml bottles	 X 16	 X 24	 X 36	 X 46
100 ml bottles	 X 8	 X 12	 X 18	 X 23
250 ml bottles	 X 6	 X 8	 X 12	 X 15
500 ml bottles	 X 3	 X 4	 X 6	 X 8
Number of trays (included)	 X 3	 X 3	 X 3	 X 3
Width of a tray	12 cm	14 cm	17 cm	17 cm
Length of a tray	30 cm	34 cm	33 cm	43 cm
Number of bottle trays (optional)	 X 1	 X 1	 X 1	 X 1
Width of bottle tray	8 cm	10 cm	12 cm	12 cm
Length of bottle tray	30 cm	34 cm	33 cm	43 cm



CertoClav Vacuum Pro Series

THE ALL-IN-ONE AUTOCLAVE

TECHNICAL DATA



Chamber volume	8 liter	12 liter	18 liter	22 liter
Product category	Laboratory horizontal table-top autoclave	Laboratory horizontal table-top autoclave	Laboratory horizontal table-top autoclave	Laboratory horizontal table-top autoclave
Order number for autoclave with integrated water tanks (stand-alone)	8500486	8500487	8500488	8500489
Order number for autoclave with fixed water connection. Including RO-System for automatic refill	8500486W	8500487W	8500488W	8500489W
Vacuuming	✓	✓	✓	✓
Drying	✓	✓	✓	✓
Liquid media sensor	✓	✓	✓	✓
User-defined program	✓	✓	✓	✓
Printer port	✓	✓	✓	✓
Ethernet	✓	✓	✓	✓
WiFi	✓	✓	✓	✓
CertoCloud Software	✓	✓	✓	✓
Product size	56 x 44 x 36 cm	61 x 50 x 43 cm	61 x 50 x 43 cm	61 x 50 x 43 cm
Diameter of chamber	17 cm	20 cm	25 cm	25 cm
Depth of chamber	32 cm	36 cm	35 cm	45 cm
Net weight	35 kg	51 kg	53 kg	55 kg
Gross weight	42 kg	59 kg	61 kg	63 kg
Carton size	65 x 45 x 42 cm	73 x 57 x 51 cm	73 x 57 x 51 cm	73 x 57 x 51 cm
Scope of delivery	1 autoclave, 5 trays, 1 tray holder, 1 tray handle, 1 power cable, 2 quick-connect drain hoses, 1 CertoCloud license	1 autoclave, 5 trays, 1 tray holder, 1 tray handle, 1 power cable, 2 quick-connect drain hoses, 1 CertoCloud license	1 autoclave, 5 trays, 1 tray holder, 1 tray handle, 1 power cable, 2 quick-connect drain hoses, 1 CertoCloud license	1 autoclave, 5 trays, 1 tray holder, 1 tray handle, 1 power cable, 2 quick-connect drain hoses, 1 CertoCloud license
Sterilizable materials	Solids, hollow materials, porous materials, liquids, prions, textiles, plastics	Solids, hollow materials, porous materials, liquids, prions, textiles, plastics	Solids, hollow materials, porous materials, liquids, prions, textiles, plastics	Solids, hollow materials, porous materials, liquids, prions, textiles, plastics
Locking mechanism	Electric door lock	Electric door lock	Electric door lock	Electric door lock
Cooling system	Air cooling	Air cooling	Air cooling	Air cooling
Screen	7" Touch screen	7" Touch screen	7" Touch screen	7" Touch screen
Chamber material	Stainless steel 304	Stainless steel 304	Stainless steel 304	Stainless steel 304
Protocol export	USB, SD-card, cloud	USB, SD-card, Cloud	USB, SD-card, cloud	USB, SD-card, cloud
Languages	English, German	English, German	English, German	English, German
Sterilization temperature	121°C to 134°C	121°C to 134°C	121°C to 134°C	121°C to 134°C
Maximum working pressure	230 kPa	230 kPa	230 kPa	230 kPa
Power supply	220V, 50/60Hz, 2400W	220V, 50/60Hz, 2400W	220V, 50/60Hz, 2400W	220V, 50/60Hz, 2400W
Certification	CE	CE	CE	CE
Warranty	2 years	2 years	2 years	2 years



CertoClav Sterilizer GmbH
Georg Grininger-Str. 37
A-4050 Traun



Phone: +437229689
Email: support@certoclav.com
Website: www.certoclav.com



Orders, dealer verification
Custom manufacture
Operating instructions and operation videos
Technical support